

Used Products

TOOTH - Crown - Metal - Retentive preparation - Super- and equigingival - Non-visible margin - Vivaglass CEM

VivaglassCEM PL

VivaglassCEM PL is a highly translucent self-curing, radiopaque glass-ionomer cement



Proxyl fluoride-free

Prophy paste without fluoride



OptraStick

Application instrument that features a flexible adhesive tip



OptraGate

Allows lips and cheeks to be retracted completely and ensures relative isolation



OptraPol

OptraPol is excellently suitable for finishing and polishing all popular composite materials in a single step



Fluor Protector

Fluor Protector is a protective fluoride varnish for desensitization and caries prophylaxis



Flowchart Vivaglass CEM

TOOTH - Crown - Metal - Retentive preparation - Super- and equigingival - Non-visible margin - Vivaglass CEM

1 The temporary is removed



The temporary is removed. If necessary, any leftover temporary cement is removed from the preparation with a polishing brush and cleaning paste free of oil and fluoride (e.g. **Proxyt fluoride-free**). Subsequently, the preparation is dried with moisture-free and oil-free air.

2 The restoration is tried in



The permanent restoration is tried in. At this stage, the shade, accuracy of fit and occlusion of the restoration are checked.

3 The restoration is pretreated



The inner surfaces of the restoration are sandblasted as directed by the manufacturer of the restorative materials.

4 The preparation is isolated and cleaned



The preparation is cleaned with a polishing brush and moisture-free and fluoride-free cleaning paste (e.g. **Proxyl fluoride-free**). Then it is rinsed with water spray. Subsequently, it is dried with air free of oil and moisture. Overdrying must be avoided.

5 Vivaglass CEM is mixed and applied



The **Vivaglass CEM** powder and liquid are mixed in a 1:1 ratio.



The luting material is applied to the restoration with a spatula or brush.

6 The restoration is seated and excess cement is removed




The restoration is seated and held in place using light constant pressure.



Once the cement has completely set, excess is removed e.g. with a scaler. The setting time is 4-6 minutes.

7 The completed restoration is finished

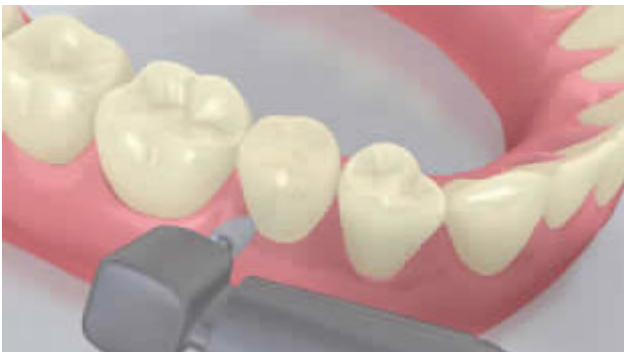

Proximal areas are adjusted with finishers and polishers. The occlusion and functional movements are checked and adjusted if necessary. The restoration margins are polished with polishers (**OptraPol**) or discs.

8 The teeth are fluoridated



A thin film of **Fluor Protector** is applied with a Vivabrush or brush and distributed evenly. The varnish is dried with an air syringe.